

TCLab User Guide

February 9, 2018

Contents

1	TCLab Overview	3
1.1	TCLab Architecture	3
1.2	Getting Started	4
1.2.1	Installation	4
1.2.2	Hardware Setup	4
1.2.3	Checking Everything Works	4
1.3	Next Steps	5
1.3.1	Course Web Sites	5
2	Accessing the Temperature Control Laboratory	6
2.1	Importing tclab	6
2.2	Using TCLab with Python's with statement	6
2.3	Reading Temperatures	7
2.4	Setting Heaters	7
2.5	Setting Maximum Heater Power	9
2.6	tclab Sampling Speed	9
2.6.1	Temperature Sampling Speed	9
2.6.2	Heater Sampling Speed	10
3	Synchronizing with Real Time	11
3.1	Simple use of tclab.clock()	11
3.2	Optional Parameters	11
3.2.1	tstep: Clock time step	11
3.2.2	tol: Clock tolerance	12
3.2.3	Fixing run time errors due to loss of clock synchronization.	12
3.3	Using tclab.clock() with TCLab	12
4	Simulation of TCLab for Offline Use	14
4.0.1	Speedup Factor with setup()	14
5	TCLab Historian	16
5.1	Basic logging	16
5.2	Accessing the Data Log from the Historian	17
5.3	Accessing log data using Pandas	18
5.4	Specifying Sources for tclab.Historian	19
6	TCLab Plotter	21
6.1	Graphics Testing	21
7	NotebookUI	23
7.0.1	Accessing past sessions	23
8	Experiment Class	24

9	Interactive and Non-blocking Operation	27
9.1	Experiments in Non-blocking Operation with <code>threading</code> Library	27
9.2	Run Class	30
9.3	Working with <code>Asyncio</code>	31
9.4	Working with <code>Tornado</code>	31
9.4.1	Adding Widgets	32

Chapter 1

TCLab Overview

The `tclab` package provides a set of Python tools for interfacing with the [BYU Temperature Control Laboratory](#). The Temperature Control Laboratory consists of two heaters and two temperature sensors mounted on an Arduino microcontroller board. Together, the `tclab` package and the Temperature Control Laboratory provide a low-cost experimental platform for implementing algorithms commonly used for process control.

1.1 TCLab Architecture


The `tclab` package is intended to be used as a teaching tool. The package provides high-level access to sensors, heaters, a pseudo-realtime clock. The package includes the following Python classes and functions:

- `TCLab()` providing access to the Temperature Control Laboratory hardware.
- `TCLabModel()` providing access to a simulation of the Temperature Control Laboratory hardware.
- `clock` for synchronizing with a real time clock.
- `Historian` for data logging.
- `Plotter` for realtime plotting.

Using these Python tools, students can create Jupyter notebooks and python codes covering a wide range of topics in process control.

- **`tclab.py`:** A Python package providing high-level access to sensors, heaters, a pseudo-realtime clock. The package includes `TCLab()` providing access to the device, `clock` for synchronizing with a real time clock, `Historian` for data logging, and `Plotter` for realtime plotting.
- **`TCLab-sketch.ino`:** Firmware for the intrinsically safe operation of the Arduino board and shield. The sketch is available at <https://github.com/jckantor/TCLab-sketch>.
- **Arduino:** Hardware platform for the Temperature Control Laboratory. TCLab is compatible with Arduino Uno, Arduino Leonardo, and compatible clones.



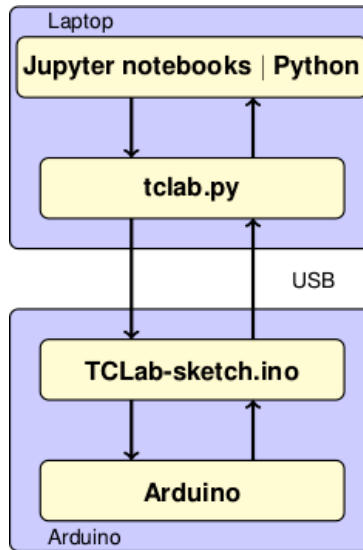

1.2 Getting Started

1.2.1 Installation

Install using

```
pip install tclab
```

To upgrade an existing installation, use the command

```
pip install tclab --upgrade
```

1.2.2 Hardware Setup

1. Plug a compatible Arduino device (UNO, Leonardo, NHduino) with the lab attached into your computer via the USB connection. Plug the DC power adapter into the wall.
2. (optional) Install Arduino Drivers.

If you are using Windows 10, the Arduino board should connect without additional drivers required.

Mac OS X users may need to install a serial driver. For Arduino clones using the CH340G, CH34G or CH34X chipset, a suitable driver can be found at <https://github.com/MPParsley/ch340g-ch34g-ch34x-mac-os-x-driver> or <https://github.com/adrianmihalko/ch340g-ch34g-ch34x-mac-os-x-driver>.

3. (optional) Install Arduino Firmware;

TCLab requires the one-time installation of custom firmware on an Arduino device. If it hasn't been pre-installed, the necessary firmware and instructions are available from the [TCLab-Sketch repository](#).

1.2.3 Checking Everything Works

Execute the following code

```
import tclab
tclab.TCLab().T1
```

If everything has worked, you should see the following output message

```
Connecting to TCLab
TCLab Firmware Version 1.2.1 on NHduino connected to port XXXX
21.54
```

The number returned is the temperature of sensor T1 in $^{\circ}\text{C}$.

1.3 Next Steps

The notebook directory provides examples on how to use the TCLab module.

1.3.1 Course Web Sites

More information, instructional videos, and Jupyter notebook examples are available at the following course websites.

- Arduino temperature control lab page <http://apmonitor.com/pdc/index.php/Main/ArduinoTemperatureControl> on the BYU Process Dynamics and Control course website.
- CBE 30338 <<http://jckantor.github.io/CBE30338/>> for the Notre Dame Chemical Process Control course website.

Chapter 2

Accessing the Temperature Control Laboratory

2.1 Importing tclab

Once installed the package can be imported into Python and an instance created with the Python statements

```
import tclab
lab = tclab.TCLab()
```

`tclab.setup()` creates an instance of a class that provides access to temperature measurements, heaters, and led on board the Temperature Control Laboratory. When called with no arguments, attempts to find a device connected to a serial port and returns a connection. An error is generated if no device is found. The connection should be closed with

```
lab.close()
```

when no longer in use. The following cell demonstrates this process, and uses the `tclab.LED()` function to flash the LED on the Temperature Control Lab for a period of 10 seconds at a 100% brightness level.

```
In [1]: import tclab
```

```
lab = tclab.TCLab()
lab.LED(100)
lab.close()
```

```
Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.
TCLab Firmware 1.3.0 Arduino Leonardo/Micro.
TCLab disconnected successfully.
```

2.2 Using TCLab with Python's with statement

The Python `with` statement provides a simple means of setting up and closing a connection to the Temperature Control Laboratory. The `with` statement establishes a context where a `tclab` instance is created, assigned to a variable, and automatically closed upon completion.

```
In [2]: import tclab
```

```
with tclab.TCLab() as lab:
    lab.LED(100)
```

```
Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.  
TCLab Firmware 1.3.0 Arduino Leonardo/Micro.  
TCLab disconnected successfully.
```

The `with` statement is likely to be the most common way to connect the Temperature Control Laboratory for most uses.

2.3 Reading Temperatures

Once a `tclab` instance is created and connected to a device, the temperature sensors on the temperature control lab can be accessed with the attributes `.T1` and `.T2`. For example, given an instance `lab`, the temperatures are accessed as

```
T1 = lab.T1  
T2 = lab.T2
```

Note that `lab.T1` and `lab.T2` are read-only properties. Any attempt to set them to a value will return a Python error.

```
In [3]: import tclab  
  
with tclab.TCLab() as lab:  
    print("Temperature 1: {0:0.2f} °C".format(lab.T1))  
    print("Temperature 2: {0:0.2f} °C".format(lab.T2))
```

```
Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.  
TCLab Firmware 1.3.0 Arduino Leonardo/Micro.  
Temperature 1: 19.62 °C  
Temperature 2: 18.65 °C  
TCLab disconnected successfully.
```

2.4 Setting Heaters

The heaters are controlled by functions `.Q1()` and `.Q2()` of a `tclab` instance. For example, both heaters can be set to 100% power with the functions

```
lab.Q1(100)  
lab.Q2(100)
```

The device firmware limits the heaters to a range of 0 to 100%. The current value of attributes may be accessed via

```
Q1 = lab.Q1()  
Q2 = lab.Q2()
```

Note that the retrieved values may be different due to the range-limiting enforced by the device firmware.

Alternatively, the heaters can also be specified with the properties `.U1` and `.U2`. Thus setting

```
lab.U1 = 100  
lab.U2 = 100
```

would set both heaters to 100% power. The current value of the heaters can be accessed as


```
print("Current setting of heater 1 is", lab.U1, "%")
print("Current setting of heater 2 is", lab.U2, "%")
```

The choice to use a function (i.e, .Q1() and .Q2()) or a property (i.e, .U1 or .U2) to set and access heater settings is a matter of user preference.

```
In [4]: import tclab
import time

with tclab.TCLab() as lab:
    print("\nStarting Temperature 1: {0:0.2f} °C".format(lab.T1),flush=True)
    print("Starting Temperature 2: {0:0.2f} °C".format(lab.T2),flush=True)

    lab.Q1(100)
    lab.Q2(100)
    print("\nSet Heater 1:", lab.Q1(), "%",flush=True)
    print("Set Heater 2:", lab.Q2(), "%",flush=True)

    t_heat = 30
    print("\nHeat for", t_heat, "seconds")
    time.sleep(t_heat)

    print("\nTurn Heaters Off")
    lab.Q1(0)
    lab.Q2(0)
    print("\nSet Heater 1:", lab.Q1(), "%",flush=True)
    print("Set Heater 2:", lab.Q2(), "%",flush=True)

    print("\nFinal Temperature 1: {0:0.2f} °C".format(lab.T1))
    print("Final Temperature 2: {0:0.2f} °C".format(lab.T2))
```

```
Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.
TCLab Firmware 1.3.0 Arduino Leonardo/Micro.
```

```
Starting Temperature 1: 19.62 °C
Starting Temperature 2: 18.65 °C
```

```
Set Heater 1: 100.0 %
Set Heater 2: 100.0 %
```

```
Heat for 30 seconds
```

```
Turn Heaters Off
```

```
Set Heater 1: 0.0 %
Set Heater 2: 0.0 %
```

```
Final Temperature 1: 21.87 °C
Final Temperature 2: 23.16 °C
TCLab disconnected successfully.
```

2.5 Setting Maximum Heater Power

The control inputs to the heaters power is normally set with functions `.Q1()` and `.Q2()` (or properties `.U1` and `.U2`) specifying a value in a range from 0 to 100% of maximum heater power.

The values of maximum heater power are specified in firmware with values in the range from 0 to 255. The default values are 200 for heater 1 and 100 for heater 2. The maximum heater power can be retrieved and set by properties `P1` and `P2`. The following code, for example, sets both heaters to a maximum power of 100.

```
In [5]: import tclab

with tclab.TCLab() as lab:
    print("Maximum power of heater 1 = ", lab.P1)
    print("Maximum power of heater 2 = ", lab.P2)

    print("Adjusting the maximum power of heater 1.")
    lab.P1 = 100

    print("Maximum power of heater 1 = ", lab.P1)
    print("Maximum power of heater 2 = ", lab.P2)
```

```
Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.
TCLab Firmware 1.3.0 Arduino Leonardo/Micro.
Maximum power of heater 1 = 200.0
Maximum power of heater 2 = 100.0
Adjusting the maximum power of heater 1.
Maximum power of heater 1 = 100.0
Maximum power of heater 2 = 100.0
TCLab disconnected successfully.
```

The actual power supplied to the heaters is a function of the power supply voltage applied to the Temperature Control Lab shield,

The maximum power applied to the heaters is a product of the settings (`P1,P2`) and of the power supply used with the TCLab hardware. The TCLab hardware is normally used with a 5 watt USB power supply capable of supply up to 1 amp at 5 volts.

The TCLab hardware actually draws more than 1 amp when both `P1` and `P2` are set to 255 and `Q1` and `Q2` are at 100%. This situation will overload the power supply and result in the power supply shutting down. Normally the power supply will reset itself after unplugging from the power mains.

Experience with the device shows keeping the sum `P1` and `P2` to a value less than 300 will avoid problems with the 5 watt power supply. If you have access to larger power supplies, then you can adjust `P1` and `P2` accordingly to achieve a wider range of temperatures.

2.6 tclab Sampling Speed

2.6.1 Temperature Sampling Speed

```
In [6]: import time
import tclab

lab = tclab.setup(connected=True)

N = 100
meas = []
```

```

with lab() as a:
    tic = time.time()
    for k in range(0,N):
        meas.append(a.T1)
    toc = time.time()

print('Reading temperature at', round(N/(toc-tic),1), 'samples per second.')

```

Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.
 TCLab Firmware 1.3.0 Arduino Leonardo/Micro.
 TCLab disconnected successfully.
 Reading temperature at 12.3 samples per second.

2.6.2 Heater Sampling Speed

```

In [7]: import time
import tclab

lab = tclab.setup(connected=True)

N = 100
with lab() as a:
    tic = time.time()
    for k in range(0,N):
        a.Q1(100)
    toc = time.time()

print('Setting heater at', round(N/(toc-tic),1), 'samples per second.')

```

Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.
 TCLab Firmware 1.3.0 Arduino Leonardo/Micro.
 TCLab disconnected successfully.
 Setting heater at 8.2 samples per second.

Chapter 3

Synchronizing with Real Time

3.1 Simple use of `tclab.clock()`

The `tclab` module includes a function `clock` for synchronizing calculations with real time. `clock(tperiod)` is an iterator that generates a sequence of equally spaced time steps from zero to `tperiod` separated by one second intervals. For each step `clock` returns time since start rounded to the nearest 10th of a second.

```
In [1]: import tclab
```

```
    tperiod = 5
    for t in tclab.clock(tperiod):
        print(t, "sec.")
```

```
0 sec.
1.0 sec.
2.0 sec.
3.0 sec.
4.0 sec.
5.0 sec.
```

`tclab.clock()` is implemented as a Python generator using the Python library `time.sleep()` function. A consequence of this implementation is that `tclab.clock()` is 'blocking' which limits its use for creating interactive demonstrations. See later sections of this user's guide for non-blocking alternatives that can be used for interactive demonstrations or GUI's.

3.2 Optional Parameters

3.2.1 `tstep`: Clock time step

An optional parameter `tstep` specifies a time step different from one second.

```
In [2]: import tclab
```

```
    tperiod = 5
    tstep = 2.5
    for t in tclab.clock(tperiod,tstep):
        print(t, "sec.")
```

```
0 sec.
2.5 sec.
```

5.0 sec.

3.2.2 tol: Clock tolerance

There are some considerations when using `clock`. First, by its nature Python is not a real-time environment. `clock` makes a best effort to stay in sync with the wall clock but there can be no guarantees. The default behavior of `clock` is to maintain long-term synchronization with the real time clock. A `RuntimeError` is raised if the difference between clock time and real time is greater than the optional parameter `tol` (default value of 0.25).

The following cell demonstrates the use of `tol` to avoid run time errors due to an intermittent calculation that exceeds the time step specified by `tstep`. In this instance, an extra `sleep` timeout of 0.2 seconds occurs at `t=1` and `t=2`. Together these cause the clock to be delayed a total of 0.4 seconds relative to real time. The default tolerance of 0.25 second would cause a run time error to be generated. Setting `tol=0.5` avoids raising that exception.

```
In [3]: import tclab
import time

tfinal = 5
tstep = 1

for t in tclab.clock(tfinal, tstep, tol=0.5):
    print(t, "sec.")
    if 0.5 < t < 2.5:
        time.sleep(1.2)
```

0 sec.
1.0 sec.
2.2 sec.
3.4 sec.
4.0 sec.
5.0 sec.

3.2.3 Fixing run time errors due to loss of clock synchronization.

If you encounter a `RuntimeError` due to loss of synchronization, there are normally two remedies depending on the underlying cause:

1. Synchronization can be lost if you attempting to multitask (say do some web browsing) while an experiment is underway. In this case you should increase the `tol` to a larger value.
2. Synchronization can be lost if the loop includes lengthy calculations, such as a computationally intensive control calculation. In the case it may be necessary to increase the `tstep` parameter.

3.3 Using `tclab.clock()` with TCLab

The `tclab.clock()` generator can be used to implement and test control and estimation algorithms using the the Temperature Control Laboratory.

```
In [4]: import tclab

tfinal = 20
tstep = 2
```

```

with tclab.TCLab() as lab:
    lab.Q1(100)
    lab.Q2(100)
    print("\nSet Heater 1 to {0:f} %".format(lab.Q1()))
    print("Set Heater 2 to {0:f} %".format(lab.Q2()))

    sfmt = "    {0:5.1f} sec:   T1 = {1:0.1f} řC    T2 = {2:0.1f} řC"
    for t in tclab.clock(tfinal, tstep):
        print(sfmt.format(t, lab.T1, lab.T2), flush=True)

```

Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.
 TCLab Firmware 1.3.0 Arduino Leonardo/Micro.

```

Set Heater 1 to 100.000000 %
Set Heater 2 to 100.000000 %
    0.0 sec:   T1 = 23.2 řC    T2 = 23.8 řC
    2.0 sec:   T1 = 23.2 řC    T2 = 23.8 řC
    4.0 sec:   T1 = 23.2 řC    T2 = 23.8 řC
    6.0 sec:   T1 = 23.2 řC    T2 = 23.8 řC
    8.0 sec:   T1 = 23.2 řC    T2 = 24.1 řC
   10.0 sec:   T1 = 23.5 řC    T2 = 23.2 řC
   12.0 sec:   T1 = 23.5 řC    T2 = 23.5 řC
   14.0 sec:   T1 = 23.5 řC    T2 = 23.8 řC
   16.0 sec:   T1 = 23.8 řC    T2 = 25.1 řC
   18.0 sec:   T1 = 23.8 řC    T2 = 25.4 řC
   20.0 sec:   T1 = 24.1 řC    T2 = 25.7 řC
TCLab disconnected successfully.

```

Chapter 4

Simulation of TCLab for Offline Use

TCLabModel replaces TCLab for occasions where the TCLab hardware might not be available. To use, include the import

```
from tclab import TCLabModel as TCLab
```

The rest of your code will work without change. Be advised the underlying model used to approximate the behavior of the Temperature Control Laboratory is an approximation to the dynamics of the actual hardware.

```
In [ ]: from tclab import TCLabModel as TCLab

        with TCLab() as a:
            print("Temperature 1: {0:0.2f} °C".format(a.T1))
            print("Temperature 2: {0:0.2f} °C".format(a.T2))
```

As an additional example.

```
In [ ]: %matplotlib inline
        from tclab import TCLabModel as TCLab
        from tclab import clock, Historian, Plotter

        with TCLab() as a:
            h = Historian(a.sources)
            p = Plotter(h, twindow=200)
            for t in clock(200):
                a.Q1(100 if t < 100 else 0)
                p.update(t)
```

4.0.1 Speedup Factor with setup()

The setup function provides control over the use of the TCLab hardware or model. If using a model, an option to run at a multiple of real-time.

A speedup of 10 or greater causes the simulation to run as fast as possible.

```
In [ ]: %matplotlib inline
        from tclab import clock, Historian, Plotter, setup

        TCLab = setup(connected=False, speedup=20)

        with TCLab() as a:
            h = Historian(a.sources)
```

```
p = Plotter(h, twindow=200)
for t in clock(200):
    a.Q1(100 if t < 100 else 0)
    h.update(t)
p.update()
```


Chapter 5

TCLab Historian

5.1 Basic logging

The `tclab.Historian` class provides data logging. Given an instance of a `TCLab` object, an historian is created with the commands

```
import tclab
lab = tclab.TCLab()
h = tclab.Historian(lab.sources)
```

The historian initializes a data log. The sources for the data log are specified in the argument to `tclab.Historian`. A default set of sources for an instance `lab` is given by `lab.sources`. (The specification for sources is described in a later section.)

The data log is updated by issuing a command

```
h.update(t)
```

Where `t` is the current clock time. If `t` is omitted the historian will calculate its own time.

```
In [1]: import tclab

        with tclab.TCLab() as lab:
            h = tclab.Historian(lab.sources)
            for t in tclab.clock(20):
                lab.Q1(100 if t <= 10 else 0)
                print("Time:", t, 'seconds')
                h.update(t)
```

```
Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.
TCLab Firmware 1.3.0 Arduino Leonardo/Micro.
```

```
Time: 0 seconds
Time: 1.0 seconds
Time: 2.0 seconds
Time: 3.0 seconds
Time: 4.0 seconds
Time: 5.0 seconds
Time: 6.0 seconds
Time: 7.0 seconds
Time: 8.0 seconds
Time: 9.0 seconds
Time: 10.0 seconds
```

```
Time: 11.0 seconds
Time: 12.0 seconds
Time: 13.0 seconds
Time: 14.0 seconds
Time: 15.0 seconds
Time: 16.0 seconds
Time: 17.0 seconds
Time: 18.0 seconds
Time: 19.0 seconds
Time: 20.0 seconds
TCLab disconnected successfully.
```

5.2 Accessing the Data Log from the Historian

Historian maintains a data log that is updated on each encounter of the `.update()` function. The list of variables logged by an Historian is given by

```
In [2]: h.columns
```

```
Out[2]: ['Time', 'T1', 'T2', 'Q1', 'Q2']
```

Individual time series are available as elements of `Historian.fields`. For the default set of sources, the time series can be obtained as

```
t, T1, T2, Q1, Q2 = h.fields
```

For example, here's how to plot the history of temperature T1 versus time from the example above.

```
In [3]: %matplotlib notebook
import matplotlib.pyplot as plt

t, T1, T2, Q1, Q2 = h.fields
plt.plot(t, T1)
plt.xlabel('Time / seconds')
plt.ylabel('Temperature / °C')
plt.grid()
```

```
<IPython.core.display.Javascript object>
```

```
<IPython.core.display.HTML object>
```

A sample code demonstrating how to plot the historian log.

```
In [4]: def plotlog(historian):
line_options = {'lw': 2, 'alpha': 0.8}
fig = plt.figure(figsize=(6, 5))
nplots = len(h.columns) - 1
t = historian.fields[0]
for n in range(1, nplots+1):
    plt.subplot(nplots,1,n)
    y = historian.fields[n]
    plt.step(t, y, where='post', **line_options)
```

```

        plt.grid()
        plt.xlabel('Time / Seconds')
        plt.ylabel(historian.columns[n])
        plt.tight_layout()

    plotlog(h)
<IPython.core.display.Javascript object>

<IPython.core.display.HTML object>

```

5.3 Accessing log data using Pandas

Pandas is a widely use Python library for manipulation and analysis of data sets. Here we show how to access the tclab.Historian log using Pandas.

The entire data history is available from the historian as the attribute `.log`. Here we show the first three rows from the log:

```

In [5]: h.log[:3]
Out[5]: [(0, 25.1, 27.03, 100.0, 0.0),
         (1.0, 25.1, 27.03, 100.0, 0.0),
         (2.0, 25.1, 27.03, 100.0, 0.0)]

```

The log can be converted to a Pandas dataframe.

```

In [6]: import pandas as pd

        df = pd.DataFrame.from_records(h.log, columns=h.columns, index='Time')
        df.head()
Out[6]:

```

	T1	T2	Q1	Q2
Time				
0.0	25.1	27.03	100.0	0.0
1.0	25.1	27.03	100.0	0.0
2.0	25.1	27.03	100.0	0.0
3.0	25.1	27.03	100.0	0.0
4.0	25.1	27.03	100.0	0.0

The following cells provide examples of plots that can be constructed once the data log has been converted to a pandas dataframe.

```

In [7]: df.plot()
<IPython.core.display.Javascript object>

<IPython.core.display.HTML object>

Out[7]: <matplotlib.axes._subplots.AxesSubplot at 0x113b1f8d0>
In [8]: df[['T1', 'T2']].plot(grid=True)
<IPython.core.display.Javascript object>

<IPython.core.display.HTML object>

Out[8]: <matplotlib.axes._subplots.AxesSubplot at 0x113f56be0>

```

5.4 Specifying Sources for tclab.Historian

An instance of `tclab.Historian` created by specifying a set of sources which are to be logged during the course of an experiment. For many cases the default sources created for an instance of `TCLab` is sufficient. However, it is possible to add additional sources which can be useful when implementing more complex algorithms for process control.

Sources are specified as a list of tuples. Each tuple has two elements. The first element is a label for the source. The second element is a function that returns a value.

The following cell shows how to create a source with the label `Power` with a value equal to the estimated heater power measured in watts. (This is created on the assumption that 100% of a maximum power of 200 corresponds to 4.2 watts).

```
In [9]: import tclab

with tclab.TCLab() as lab:
    sources = [
        ('T1', lambda: lab.T1),
        ('Power', lambda: lab.P1*lab.U1*4.2/(200*100))
    ]
    h = tclab.Historian(sources)
    for t in tclab.clock(20):
        lab.Q1(100 if t <= 10 else 0)
        print("Time:", t, 'seconds')
        h.update(t)
```

Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.

TCLab Firmware 1.3.0 Arduino Leonardo/Micro.

Time: 0 seconds

Time: 1.0 seconds

Time: 2.0 seconds

Time: 3.0 seconds

Time: 4.0 seconds

Time: 5.0 seconds

Time: 6.0 seconds

Time: 7.0 seconds

Time: 8.0 seconds

Time: 9.0 seconds

Time: 10.0 seconds

Time: 11.0 seconds

Time: 12.0 seconds

Time: 13.0 seconds

Time: 14.0 seconds

Time: 15.0 seconds

Time: 16.0 seconds

Time: 17.0 seconds

Time: 18.0 seconds

Time: 19.0 seconds

Time: 20.0 seconds

TCLab disconnected successfully.

```
In [10]: import pandas as pd
```

```
df = pd.DataFrame.from_records(h.log, columns=h.columns, index='Time')
df.head()
```

```
Out[10]:
```

	T1	Power
Time		
0.0	26.39	4.2
1.0	26.39	4.2
2.0	26.39	4.2
3.0	26.39	4.2
4.0	26.39	4.2

Chapter 6

TCLab Plotter

When operating in a Jupyter Notebook, a Plotter can be used together with the Historian.

```
h = Historian(a)
p = Plotter(h, tfinal)
```

where `a` is a TCLab instance as before and the optional parameter `tfinal` provides an initial scaling of the time axes. Each call to `p.update()` will automatically update both the historian and the plot.

```
In [1]: %matplotlib notebook
        from tclab import TCLab, clock, Historian, Plotter

        with TCLab() as a:
            h = Historian(a.sources)
            p = Plotter(h, twindow=200)
            for t in clock(200):
                a.Q1(100 if t < 100 else 0)
                p.update(t)
```

```
Arduino Leonardo connected on port /dev/cu.usbmodemWUAR1 at 115200 baud.
TCLab Firmware 1.3.0 Arduino Leonardo/Micro.
```

```
<IPython.core.display.Javascript object>
```

```
<IPython.core.display.HTML object>
```

```
TCLab disconnected successfully.
```

6.1 Graphics Testing

Notebook to test graphics during development.

```
In [2]: %matplotlib notebook
        from tclab import setup
        from tclab import Historian, Plotter, clock
        import time
```

```
tic = time.time()
lab = setup(connected=False, speedup=10)
with lab() as a:
    h = Historian(a.sources)
    p = Plotter(h, 200, layout= (('T1', 'T2'), ('Q1', 'Q2')))
    for t in clock(200):
        a.U1 = 80
        p.update(t)
toc = time.time()

print(toc-tic, 'seconds')
```

Simulated TCLab

<IPython.core.display.Javascript object>

<IPython.core.display.HTML object>

TCLab Model disconnected successfully.
12.71963119506836 seconds

Chapter 7

NotebookUI

The `tclab.gui` module supplies a graphical interface to the Temperature Control Laboratory.

```
In [1]: from tclab.gui import NotebookUI
```

```
In [2]: %matplotlib notebook
```

```
In [3]: interface = NotebookUI()
```

```
In [4]: interface.gui
```

A Jupyter Widget

7.0.1 Accessing past sessions

Once you have finished the experiment, you can see what sessions the historian stored as follows:

```
In [5]: interface.historian.get_sessions()
```

```
-----  
  
AttributeError                                Traceback (most recent call last)  
  
  <ipython-input-5-659f5f7e1ad4> in <module>()  
----> 1 interface.historian.get_sessions()  
  
AttributeError: 'NotebookUI' object has no attribute 'historian'
```

The historian can load data from one of the previous sessions. Note that this will overwrite the data currently stored in the historian.

```
In [ ]: interface.historian.load_session(1)
```

Once the data have been loaded, the historian will support all the same commands as if only one session was used.

```
In [ ]: interface.historian.log
```


Chapter 8

Experiment Class

```
In [1]: from tclab import runexperiment
```

```
In [2]: %matplotlib notebook
```

```
In [3]: def onoffcontroller(t, lab):  
        lab.Q1(0 if t < 5 else 100)
```

```
In [4]: experiment = runexperiment(onoffcontroller, connected=False, time=100, speedup=10)
```

Simulated TCLab

<IPython.core.display.Javascript object>

<IPython.core.display.HTML object>

TCLab Model disconnected successfully.

```
In [5]: experiment.historian.log
```

```
Out[5]: [(0, 21.114271687292348, 21.136557764983188, 0, 0),  
        (1, 20.989368089438035, 21.002784474310026, 0, 0),  
        (2, 20.995700325216735, 20.79449784532808, 0, 0),  
        (3, 20.95021966154697, 20.808941623501386, 0, 0),  
        (4, 20.930766355723883, 20.96593591197957, 0, 0),  
        (5, 20.936535428781973, 21.036142472121153, 0, 0),  
        (6, 20.91787552408186, 20.86723366607012, 100, 0),  
        (7, 20.88244225305503, 21.221179444040317, 100, 0),  
        (8, 21.421053184788743, 21.194534902582085, 100, 0),  
        (9, 21.355918262585735, 20.944441931004164, 100, 0),  
        (10, 20.939346577272488, 20.85249048840826, 100, 0),  
        (11, 21.293909172212985, 20.885712402297678, 100, 0),  
        (12, 21.418291938150755, 21.16923119425466, 100, 0),  
        (13, 21.495768879101302, 21.121956476044403, 100, 0),  
        (14, 21.91453057133991, 20.716846733269314, 100, 0),  
        (15, 21.713962135869807, 21.12718006542719, 100, 0),  
        (16, 22.484961718780234, 21.225676788916306, 100, 0),  
        (17, 22.32081303475573, 21.160418468648267, 100, 0),
```

(18, 22.869546717455897, 20.907974280417847, 100, 0),
(19, 22.931785538132353, 21.215296996953544, 100, 0),
(20, 23.208130725516792, 20.59571158741274, 100, 0),
(21, 22.965400032019943, 21.05860705963216, 100, 0),
(22, 23.459051960378872, 21.22957341575619, 100, 0),
(23, 23.22090474187922, 21.557769880711323, 100, 0),
(24, 24.189015585930978, 21.27653739408498, 100, 0),
(25, 24.012767830097026, 21.18833546964699, 100, 0),
(26, 24.3987951033398, 21.296745611212092, 100, 0),
(27, 24.996892462628992, 21.077275044590827, 100, 0),
(28, 25.05615676638348, 21.0349736469255, 100, 0),
(29, 24.853770446471643, 21.03000433705572, 100, 0),
(30, 25.785163927472738, 21.129450998361644, 100, 0),
(31, 26.333833450814193, 20.99187156503794, 100, 0),
(32, 26.16400141106771, 21.157186636725203, 100, 0),
(33, 26.523337085749663, 21.372817700106566, 100, 0),
(34, 26.80131870216762, 21.91084688249117, 100, 0),
(35, 27.22293927686771, 21.409085032771134, 100, 0),
(36, 27.37157026264135, 21.618728734635408, 100, 0),
(37, 27.53842911416146, 21.80416591846843, 100, 0),
(38, 27.760688657574104, 21.263426600560955, 100, 0),
(39, 28.337488852702453, 21.600046538474842, 100, 0),
(40, 29.03441457434998, 21.989108598841856, 100, 0),
(41, 29.031260441916885, 21.590822081318223, 100, 0),
(42, 29.606725859873553, 21.996239818767073, 100, 0),
(43, 29.789947333106888, 21.517566998525833, 100, 0),
(44, 30.222109393953957, 21.783042199378375, 100, 0),
(45, 30.73259605665824, 21.848586846428482, 100, 0),
(46, 30.522275399508363, 22.038279492463413, 100, 0),
(47, 31.05369807941251, 22.125503178438837, 100, 0),
(48, 31.45127095668413, 22.05064774445467, 100, 0),
(49, 31.542652686032106, 22.11930176964894, 100, 0),
(50, 32.15947283709359, 22.08527341932343, 100, 0),
(51, 32.57879615212834, 22.002545318785007, 100, 0),
(52, 32.63246371256205, 22.224998110048492, 100, 0),
(53, 33.27716112687957, 22.058858818011426, 100, 0),
(54, 33.30157120876311, 22.417030156650505, 100, 0),
(55, 33.206360625826846, 22.626680851638525, 100, 0),
(56, 33.72311365767237, 22.42776780703317, 100, 0),
(57, 33.94024059283053, 22.40763618975859, 100, 0),
(58, 34.19319480304237, 22.529969928862027, 100, 0),
(59, 34.762254310835296, 22.219713588090052, 100, 0),
(60, 34.964690022686874, 22.467076360258577, 100, 0),
(61, 35.36469614861533, 22.914230766735486, 100, 0),
(62, 35.57224302671094, 22.524005558951234, 100, 0),
(63, 35.77769304515034, 22.62697158082207, 100, 0),
(64, 36.598684087429, 22.639552422850475, 100, 0),
(65, 36.48756343520329, 23.0422224423655, 100, 0),
(66, 37.02834749115227, 22.823455484534264, 100, 0),
(67, 37.56328311999661, 22.915257421135657, 100, 0),
(68, 37.7990785380291, 22.86157148330309, 100, 0),
(69, 38.03301907855134, 22.755446520644753, 100, 0),
(70, 38.15830040730779, 23.060959365611083, 100, 0),
(71, 38.19112008767047, 23.51999032676287, 100, 0),

(72, 39.025828250588255, 23.092923897284628, 100, 0),
 (73, 38.932639874359914, 23.470496027626155, 100, 0),
 (74, 39.179538406104, 23.31962217322741, 100, 0),
 (75, 39.39210484339698, 23.54251946180979, 100, 0),
 (76, 39.83268776885665, 23.602422574686646, 100, 0),
 (77, 40.22172905940857, 23.250843921124137, 100, 0),
 (78, 40.14358943198019, 23.706291880666324, 100, 0),
 (79, 40.462458667404, 23.186335970556083, 100, 0),
 (80, 41.03479735771997, 23.53755270074161, 100, 0),
 (81, 41.24029241319539, 23.544017636982016, 100, 0),
 (82, 41.50890709240405, 23.88069977409174, 100, 0),
 (83, 41.923461581588576, 23.852920651090596, 100, 0),
 (84, 42.11478701750082, 23.544152638178534, 100, 0),
 (85, 42.40750886693584, 23.519830014780425, 100, 0),
 (86, 42.511198372140946, 23.337989675833647, 100, 0),
 (87, 42.41241553764143, 23.704944974921958, 100, 0),
 (88, 43.36494261271852, 24.166321401181747, 100, 0),
 (89, 43.1979805962873, 24.031787021813287, 100, 0),
 (90, 43.53156038335274, 24.11065802969447, 100, 0),
 (91, 43.85112079557571, 24.160546969157572, 100, 0),
 (92, 44.38170967732593, 24.253999556475673, 100, 0),
 (93, 44.83738785143763, 23.960311122247592, 100, 0),
 (94, 44.32127846613241, 24.30052823828633, 100, 0),
 (95, 44.916852219833075, 24.219136910475214, 100, 0),
 (96, 45.18349432758581, 24.493899341406205, 100, 0),
 (97, 45.641524150698515, 24.358971405466804, 100, 0),
 (98, 45.90942032182922, 24.48433405906429, 100, 0),
 (99, 45.82455447968922, 24.252748076206704, 100, 0),
 (100, 45.838597826584234, 24.212578698030825, 100, 0)]

Chapter 9

Interactive and Non-blocking Operation

The following sections in this notebook demonstrate methods for interacting with TCLab, for building non-blocking implementations of a control loop, and for various experiments and tests with the package.

9.1 Experiments in Non-blocking Operation with threading Library

The current implementation of

```
In [1]: def bar():
        clock.send(None)

        def clock(tperiod):
            tstart = time.time()
            tfinish = tstart + tperiod
            t = 0
            while t + tstart < tfinish:
                z = yield t
                t += 1

        def bar():
            clock.send(2)

In [2]: import threading, time
        import datetime

        next_call = time.time()
        k = 0

        def foo():
            global next_call, k
            if k < 5:
                print(k, datetime.datetime.now())
                next_call = next_call+1
                threading.Timer( next_call - time.time(), foo ).start()
                k += 1
            else:
                print(k, "Last Call")

        foo()
```

```
In [3]: from tclab import TCLabModel, Historian, Plotter
import threading, time
```


```
tstep = 1
tperiod = 20
```

```
tstart = time.time()
tfinish = tstart + tperiod
tnext = tstart
```

```
a = TCLabModel()
h = Historian(a.sources)
p = Plotter(h,20)
a.U1 = 100
```

```
def tasks(tnext):
    global tnext, tfinish, tstep
    p.update(tnext-tstart)
    tnext = tnext + tstep
    if tnext <= tfinish:
        threading.Timer(tnext-time.time(), update).start()
    else:
        a.close()
```

```
update()
```

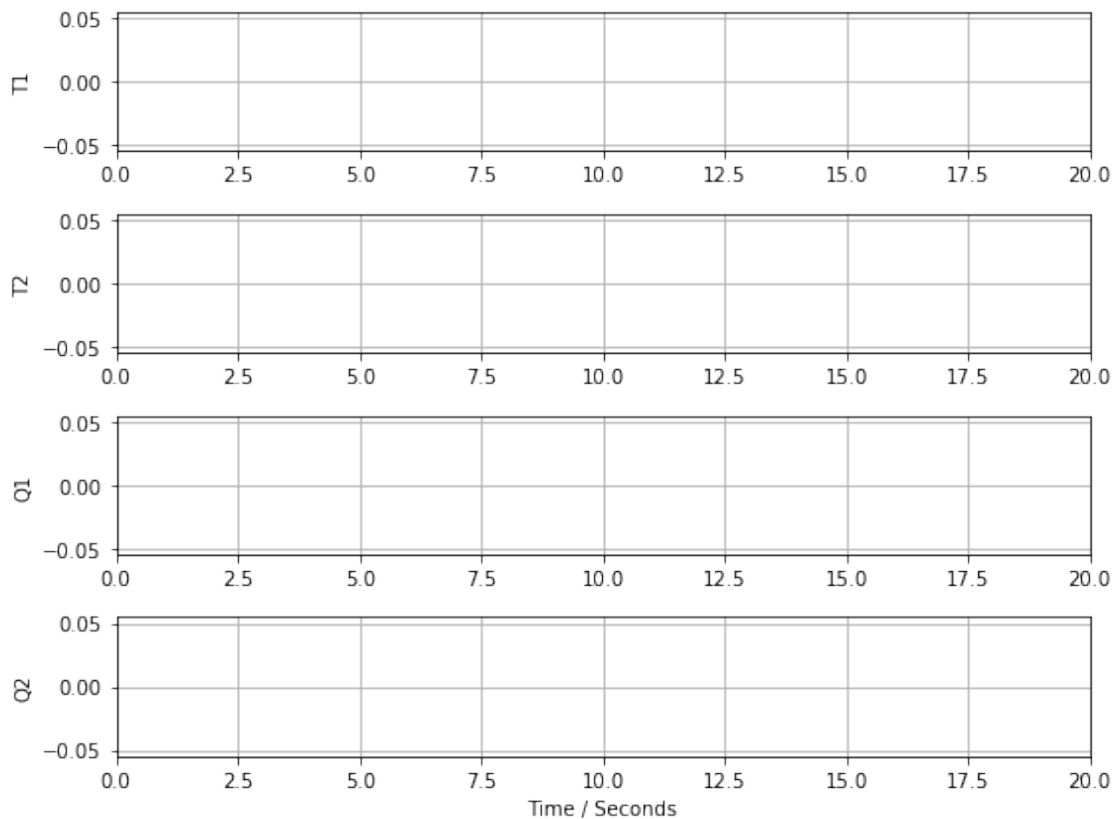

```
File "<ipython-input-3-f49db311273e>", line 18
global tnext, tfinish, tstep
^
```

SyntaxError: name 'tnext' is parameter and global

```
1 2018-02-09 20:07:40.226656
2 2018-02-09 20:07:41.221889
3 2018-02-09 20:07:42.226601
4 2018-02-09 20:07:43.224978
5 Last Call
```

```
In [ ]: %matplotlib notebook
```

```
import time
from threading import Timer
from tclab import setup, Historian, Plotter
```

```
lab = setup(connected=False, speedup=1)
a = lab()
h = Historian(a.sources)
p = Plotter(h)
```

```
SP = 40
```

```
tstart = time.time()
def loop():
    PV = a.T1
    MV = 100 if PV < SP else 0
    a.U1 = MV
    p.update(time.time()-tstart)
```

```
for t in range(0,100):
    Timer(t, loop).start()
Timer(100,a.close).start()
```

```
In [ ]: SP = 20
```

```
In [ ]: import threading, time, datetime
```

```
def loop():
    yield
    print(datetime.datetime.now())
    threading.Timer(1000, lambda: next(loop_gen)).start()
```

```
loop_gen = loop()
next(loop_gen)
```

9.2 Run Class

```
In [ ]: from threading import Timer
import time
import tclab

class Run(object):
    def __init__(self, function, tfinal, tinterval=1):
        self.lab = tclab.TCLab()
        self.tfinal = tfinal
        self.tinterval = tinterval
        self.function = function
        self._timer = None
        self.tstart = time.time()
        self.tnow = self.tstart
        self.is_running = False
        self.start()

    def _run(self):
        """Start a new timer, then run the callback."""
        self.is_running = False
        self.start()
        self.function(self.lab, self.tnow)

    def start(self):
        if not self.is_running:
            self.tnow = time.time() - self.tstart
            if self.tnow < self.tfinal:
                self._timer = Timer(self.tinterval - self.tnow % self.tinterval, self._run)
            else:
                self._timer = Timer(self.tinterval - self.tnow % self.tinterval, self.stop)
            self._timer.start()
            self.is_running = True

    def stop(self):
        print("")
        self.lab.close()
        self._timer.cancel()
        self.is_running = False

In [ ]: from ipywidgets import interact

SP = 40
Kp = 15

def loop(lab, t):
    PV = lab.T1
    MV = Kp*(SP-PV)
    lab.U1 = MV
    print("\r{0:8.2f}    {1:6.2f}    {2:6.0f}".format(t,PV,MV), end='')

expt = Run(loop, 200, 1)

In [ ]: import tclab
```

```

SP = 90

def ControlLoop(lab, t):
    PV = lab.T1
    MV = 100 if PV < SP else 0
    lab.U1 = MV
    print(round(t,4), PV, MV)
    p.update(t)

lab = tclab.TCLab()
h = tclab.Historian(lab.sources, dbfile=None)
p = tclab.Plotter(h)
expt = PeriodicCallback(lab, ControlLoop, 10, 2)

```

9.3 Working with Asyncio

```

In [ ]: import asyncio

async def slow_op(n):
    await asyncio.sleep(n)
    print("Slow Op:", n)

async def main():
    await asyncio.wait([slow_op(3),slow_op(2),slow_op(1)])

loop = asyncio.get_event_loop()
loop.run_until_complete(main())

```

9.4 Working with Tornado

This is an experiment to build a non-blocking event loop for TCLab. The main idea is to implement the main event loop as a generator, then use Tornado's non-blocking timer to send periodic messages to the generator.

```

In [ ]: %matplotlib inline
import tornado
import time
from tclab import setup, Historian, Plotter

SP = 40
Kp = 10

def update(lab):
    t = 0
    h = Historian(lab.sources)
    p = Plotter(h,120)
    while True:
        PV = lab.T1
        MV = Kp*(SP-PV)
        lab.U1 = MV
        p.update(t)
        yield

```



```

        t += 1

lab = setup(connected=True)
a = lab()
update_gen = update(a)
timer = tornado.ioloop.PeriodicCallback(lambda: next(update_gen), 1000)
timer.start()

In [ ]: timer.stop()
        a.close()

```

9.4.1 Adding Widgets

tclab.clock is based on a generator, which maintains a single thread of execution. One consequence is that there is no interaction with Jupyter widgets.

```

In [ ]: from ipywidgets import interactive
        from IPython.display import display
        from tclab import clock

Kp = interactive(lambda Kp: Kp, Kp = 12)
display(Kp)

for t in clock(10):
    print(t, Kp.result)

In [ ]: import tornado
        from ipywidgets import interactive
        from IPython.display import display
        from tclab import TCLab, Historian, Plotter

Kp = interactive(lambda Kp: Kp, Kp = (0,20))
SP = interactive(lambda SP: SP, SP = (25,55))
SP.layout.height = '500px'

def update(tperiod):
    t = 0
    with TCLab() as a:
        h = Historian(a.sources)
        p = Plotter(h)
        while t <= tperiod:
            yield
            p.update(t)
            display(Kp)
            display(SP)
            a.U1 = SP.result
            t += 1
        timer.stop()

update_gen = update(20)
timer = tornado.ioloop.PeriodicCallback(lambda: next(update_gen), 1000)
timer.start()

In [ ]: from ipywidgets import interactive
        from tclab import setup, clock, Historian, Plotter

```

```

def proportional(Kp):
    MV = 0
    while True:
        PV, SP = yield MV
        MV = Kp*(SP-PV)

def sim(Kp=1, SP=40):
    controller = proportional(Kp)
    controller.send(None)

    lab = setup(connected=False, speedup=20)
    with lab() as a:
        h = Historian(a.sources)
        p = Plotter(h,200)
        for t in clock(200):
            PV = a.T1
            MV = controller.send([PV,SP])
            a.U1 = MV
            h.update()
            p.update()

interactive_plot = interactive(sim, Kp=(0,20,1), SP=(25,60,5), continuous_update=False);
output = interactive_plot.children[-1]
output.layout.height = '500px'
interactive_plot

In [ ]: timer.stop()

```